




**THE No. 1 System
for Solvent Safety**

NEW: "Arnold" **Ball-Valve Funnel V2.0** **with Hinged Lid**

Strong Protection when collecting
Liquid Waste!

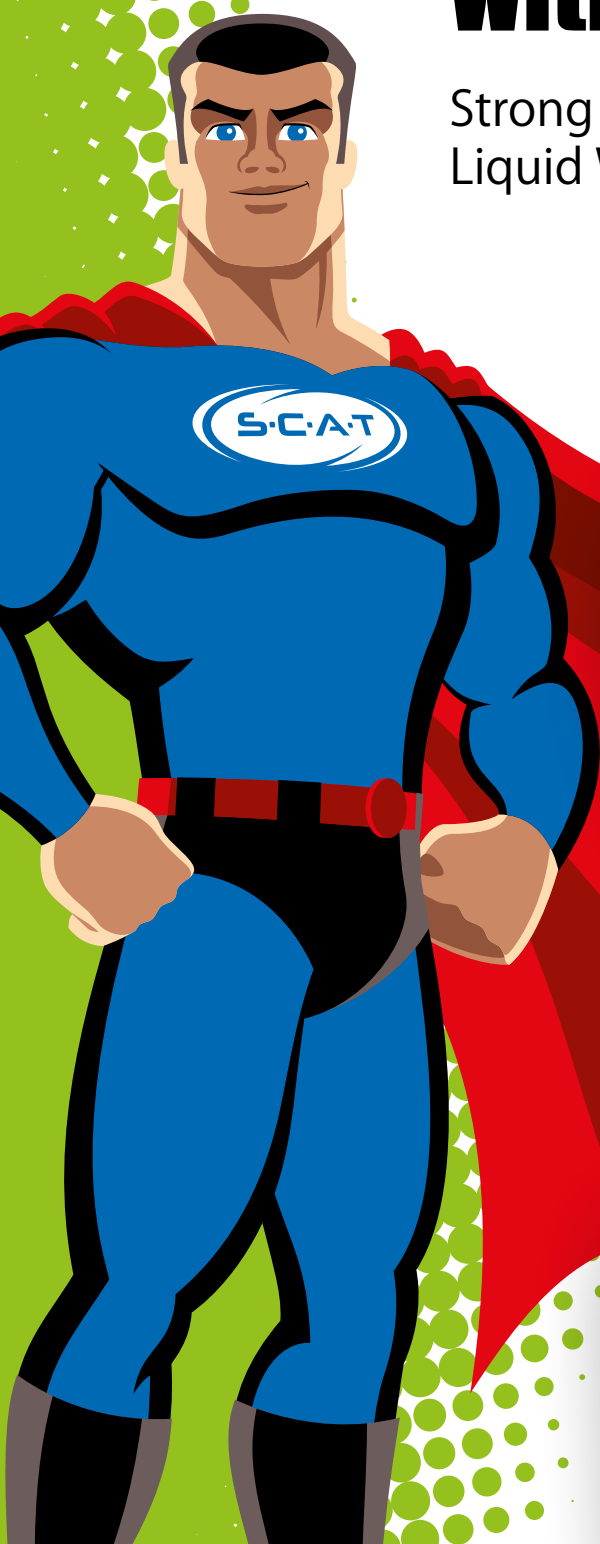

www.scateurope.com

Clever Technology, Strong Protection!

Safety Funnel "Arnold"

✓ **NEW: Improved Closing**

Comfortable handling: easy, quick opening and closing with "click"-closure.

✓ **NEW: Hinged Lid with Sealing Ring**

For secure closure, when the funnel is not in use.

✓ **NEW: Stable Hinge**

A double-screw hinge, constructed of rust-free stainless steel: for long life and stability.

✓ **NEW: Optimized Spray Protection**

Air ducts directed downwards ensure for secure ventilation and cause any spray created to be returned to the funnel.

✓ **NEW: Optimized Dirt Sieve**

For collecting mixing magnets or large dirt particles. Laboratory bottles can also be placed here to dry off. The hole size has been optimized, in order to improve the running-off of chemicals.

Both sieve and spray protection can be removed for cleaning.

✓ **Ball-Valve**

The valve only opens up during filling, and then closes automatically afterwards - for a big "Plus" regarding safety!

✓ **Various Screw-Threads**



Suitable for different sizes of collecting canister (see Table on Page 3).

✓ **Safety Lance**



Ensures for safe and good flowing of sensitive chemicals, without any "stuttering".

✓ **NEW: Robust Design**

The new design provides for a high degree of stability as well as long life, and prevents any deformation of the funnel from occurring.

Information for Ordering

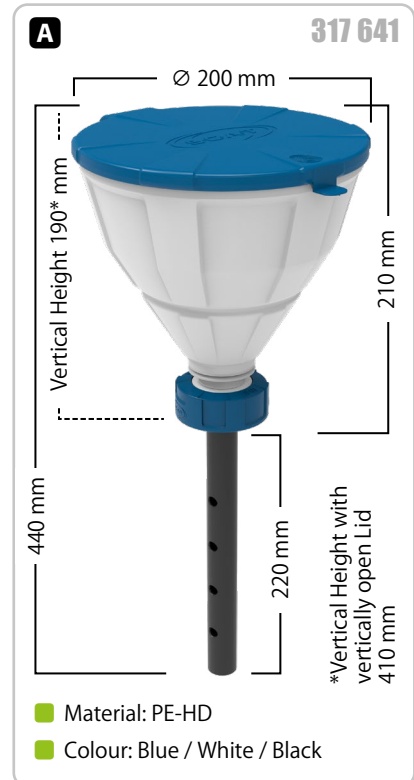


Fig.	Part No.	Screw-Thread of Container	Lance	Electrically conductive	Delivery with Adapter
-	317 642	GL 45	●	-	-
-	317 622	GL 45	●	●	-
-	317 638	GL 45	-	●	-
-	317 649	S 50	●	-	●
-	317 629	S 50	●	●	●
-	317 644	S 51	●	-	●
-	317 624	S 51	●	●	●
-	317 632	B 53	●	-	●
-	317 645	S 55	●	-	-
-	317 625	S 55	●	●	-
A	317 641	S 60 / 61	●	-	-
B	317 621	S 60 / 61	●	●	-
-	317 646	S 65	●	-	●
-	317 626	S 65	●	●	●
-	317 648	S 70 / 71	●	-	●
-	317 628	S 70 / 71	●	●	●
-	317 647	B 83	●	-	●
-	317 627	B 83	●	●	●
-	317 643	S 90	●	-	●
-	317 623	S 90	●	●	●
-	317 635	S 95	●	●	●
-	317 630	Dual Screw-Thread R2" BSP/ G2" + 2" TriSure	-	●	-
Models without ball valve, for highly viscous liquids					
-	317 633	GL 45	-	●	-
-	317 634	S 60 / 61	-	●	-



Information about "Delivery with Adapter":

The safety funnel "Arnold" is available in various container screw-thread sizes. In order that we may produce in an environmentally-friendly and economic manner, certain variants are delivered together with a screw-thread adapter. Those to which this applies can be seen in the column "Delivery with Adapter".



* The total height can vary, dependent upon the thread adapter with which the funnel has been delivered.





S.C.A.T. Europe GmbH

Waldecker Straße 7
64546 Mörfelden-Walldorf
Deutschland

 **+49 (0) 61 05 / 30 55 86 - 0**

 **+49 (0) 61 05 / 30 55 86 - 99**

 **info@scat-europe.com**

 **www.scat-europe.com**

S.C.A.T. Americas, Inc.

1704 East Blvd., Suite 101
Charlotte, NC 28203
USA

 **+1.704.438.9477**

 **+1.704.373.2603**

 **info@scat-americas.com**

 **www.scat-americas.com**



Exonder S.r.l.

 **Via Torre, 2 - 28040 Borgo Ticino (NO) - Italy**

 **+39 0321 908466**

 **customer@exonder.it**

 **https://www.exonder.it**